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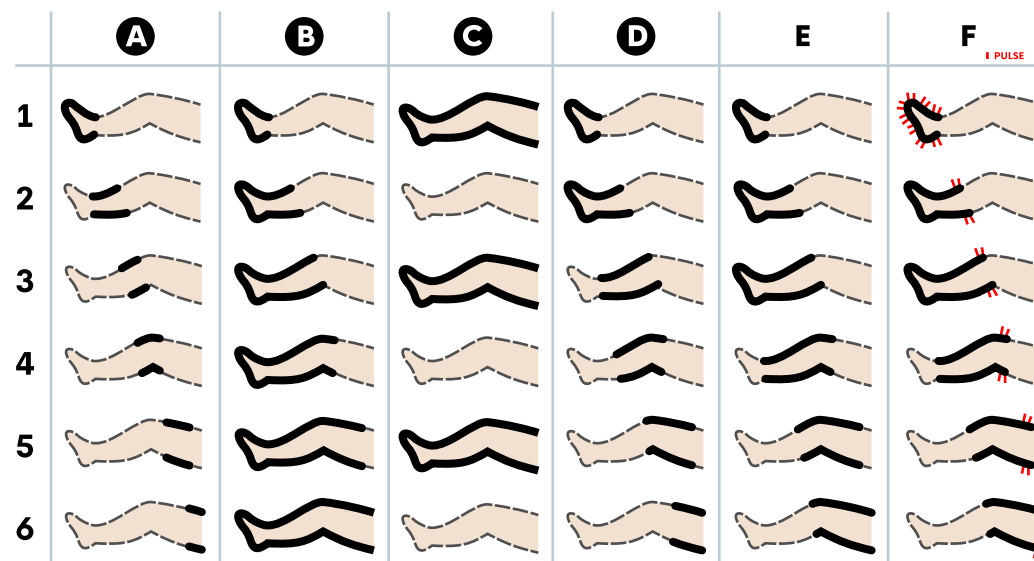
# CLEANING AND MAINTENANCE


- Turn off the compressor before cleaning
- Do not wash the compressor and cuffs
- Wipe with a dry or damp cloth.
- Keep in a clean and dry place.
- Protect the compressor and the equipment from fire, direct sunlight or heat sources.

## WARNING

- Before the first use, read the user manual included carefully
- Keep out of reach of children.
- Do not unzip the cuffs when the massage program is on.
- Protect the compressor from water and moisture.
- Before switching the compressor on, make sure that it is on a flat surface.
- Do not move the cuffs when using the device.
- The device cannot be repaired or maintained in any way during use.
- Consult your doctor or therapist before use.
- This product is not intended to diagnose, treat, or cure any diseases.
- The compressor contains a built-in Li-ion battery. Follow applicable regulations and guidelines when disposing.
- Contraindications include (but are not limited to): acute inflammation, skin disease, cardiac arrhythmia, erysipelas, deep vein thrombosis, pulmonary edema, unstable high blood pressure, patients with a pacemaker.
- Do not use with damaged liners. Replace damaged cuffs immediately.
- Use only with original Spophy Air Recovery components.
- Do not use in a humid environment.
- Do not attempt to repair or disassemble the device yourself.
- To unplug, grasp only the charging adapter.
- Overcharging the battery may shorten its service life.
- In case of damage, immediately disconnect the device from the power supply and stop using it.

spophy

TRAIN  
PHYSIO  
SLEEP  
REPEAT

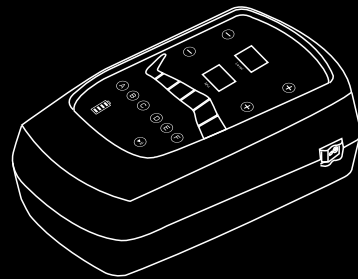Distributor: Spophy s.r.o.  
www.spophy.com  
Made in P.R.C.


spophy

# AIR RECOVERY BOOTS


TYPE: SPOBOOTS6

## INSTRUCTIONS FOR USE


- Connect according to illustration

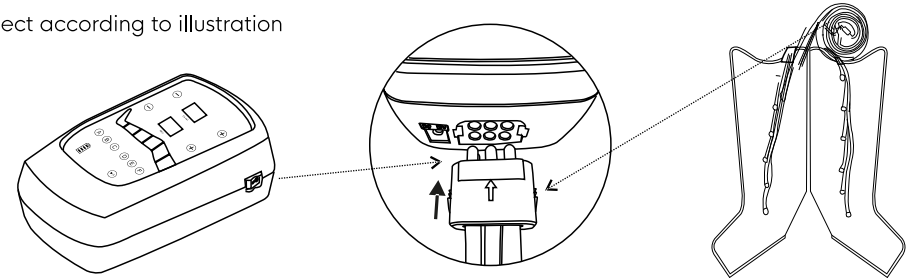- After use and when preparing of the compressor for storage, make sure that the hose connector is disconnected while pressing both safety locks. This will prevent damage to the connector plug.

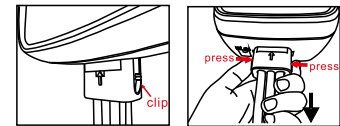- Disconnect the hoses from the cuff one at a time. Avoid disconnecting all hoses at the same time. For travel and to optimize the space in the bag, you can disconnect the compressor from the cuffs and pull out all the connectors. Frequent disconnection of the hoses from the cuffs is undesirable.

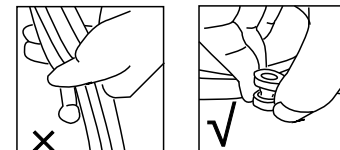- Connect the hoses to the cuffs as shown so that the longest hose is connected to the farthest segment in the ankle area.
- When disconnecting the hoses from the cuff, always grasp the plastic end of the hose and hold the lower part under the connector and, carefully pull it out, do not pull on the hose.

# COMPRESSOR SPECIFICATIONS




Dimensions: 245 mm x 150 mm x 90 mm  
 Power consumption: 35 W  
 Charger input: AC 110 V ~ 240 V; 50 Hz ~ 60 Hz  
 Charger output: DC 12.6 V; 2 A  
 Pressure range: 30 ~ 240 mmHg (± 30 mmHg)

## SETUP GUIDE

Before use:  
 Check the battery charge status by pressing the green compressor start button.  
 If the compressor does not start, it means that the battery is flat.  
 Use the charger to charge the battery.


 Only use the original charger provided in the package, other types of chargers and components may permanently damage the device.

- Securely connect the air hoses to the cuffs.
- Securely connect the main hose connector to the device.
- Put on your cuffs comfortably.
- Press the power button on the right side of the compressor - this will put the device into standby mode, from which you can further preset the desired massage. Press this button again to turn off the device.
- Set the desired pressure level using the  and  buttons.  
 Note: The default level is 8 (85 mmHg).  
 1~4 30 - 75 mmHg  
 5~8 85 - 130 mmHg  
 9~12 140 - 185 mmHg  
 13~16 195 - 240 mmHg
- Set the required treatment time using the  and  buttons.  
 Note: the time range is 5 - 99 minutes. The default time is 30 minutes. One short press increases / decreases the time by 1 minute. Hold down the button to ensure that a larger change in the desired treatment time is set faster. After the set treatment time runs out, the compressor switches off automatically.
- Select the desired massage area using the buttons showing the leg (chamber) segments. Note: the default setting is all 6 segments with the light on. Individual segments during treatment can only be switched on or off in pause mode. When the segment button light is on, it is turned on.
- Select the desired massage program using the buttons. Note: the default program is A.
- Press the start / pause button to start 
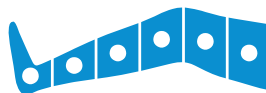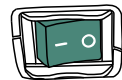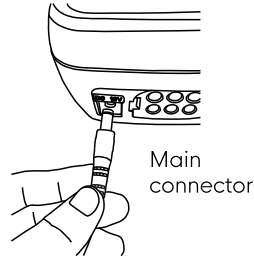
# MESSAGE PROGRAMS

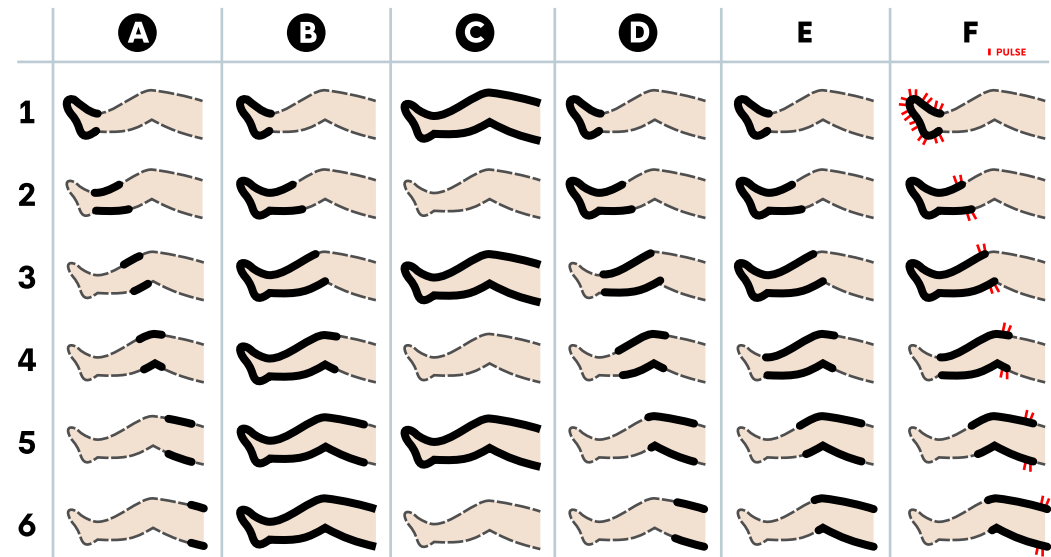- In this program, only one chamber is inflated at a time. The chamber is deflated when the required pressure is reached, then it is inflated and the next chamber is deflated. The cycle repeats.
- During this program, one chamber remains inflated and the other chambers are gradually inflated separately. After inflating all 6 segments, all chambers are deflated simultaneously. The cycle repeats.
- In this program, all chambers are inflated and deflated simultaneously. The cycle repeats.
- In this program, the first 2 chambers are gradually inflated and then the pressure of the 2 chambers is increased. In this program, the 1st chamber is inflated and holds its pressure, then the 2nd chamber is inflated and both hold their pressure, when the 3rd chamber is inflated, the 1st chamber is deflated and only the 2nd and 3rd chambers hold their pressure. In this way, the chambers are gradually inflated up to the final, 6th chamber. The cycle repeats.
- In this program, the first 3 chambers are gradually inflated and then the pressure of the 3 chambers is increased. The 1st chamber inflates and holds its pressure, the 2nd chamber inflates and the pressure is held in the 1st and 2nd chambers, the 3rd chamber inflates and the pressure is held in the 1st, 2nd and 3rd chambers. The 4th chamber is inflated while the 1st chamber is deflated and the pressure is maintained in the 2nd, 3rd and 4th chambers. In this way, the chambers are gradually inflated up to the final, 6th chamber. The cycle repeats.
- With this program, chamber 1 inflates and pulsates at the required pressure. When inflating the 2nd chamber, the 1st chamber maintains a constant pressure. The 2nd chamber starts pulsating when the required pressure is reached. When inflating the 3rd chamber, the 1st and 2nd chambers maintain a constant pressure. When the 4th chamber is inflated, the 1st chamber is deflated and the 2nd and 3rd chambers hold their pressure. When the desired pressure is reached, the 4th chamber begins to pulsate. In this way, the chambers are gradually inflated up to the final, 6th chamber. The cycle repeats. Note: pulsing means quickly inflating and deflating 3 times.

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